Barcodeleser

Anleitung

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1.7 Operation

Adjustment to Microplate Barcode Label

Power up instrument with pipetting head inserted but nothing else mounted to the pipetting head.

NOTE

Barcode reading requires any tips, liquid handling adapters, CyBi®-RoboTipTrays or cover magazine not being mounted to the head. The barcode reading library scripts will verify this and give an error if anything is mounted.

Start CyBio® Composer. Run empty method to calibrate drives. Add script "Mirror adjustment" from the library service "CyBi®-FeliX barcode reader" to a new method. Run method and select deck position for mirror adjustment. The instrument will move the barcode reader to the selected deck position and start reading repeatedly. The barcode reader's status LED (lower LED) will be continuously illuminated green, its decode LED (upper LED) will be continuously illuminated orange.

Fig. 2: Barcode reader's status LED continuously illuminated green, its decode LED continuously illuminated orange

Open blind.

CAUTION

Risk of injury to the eyes by laser light!


Do not stare into red laser beam or its reflection. The barcode reader is a laser class II device.

NOTE

Observe comprehensive information of the barcode reader manual, especially chapter "Safety Notices".

Turn mirror until the red reading line appears approximately in the middle of the barcode.

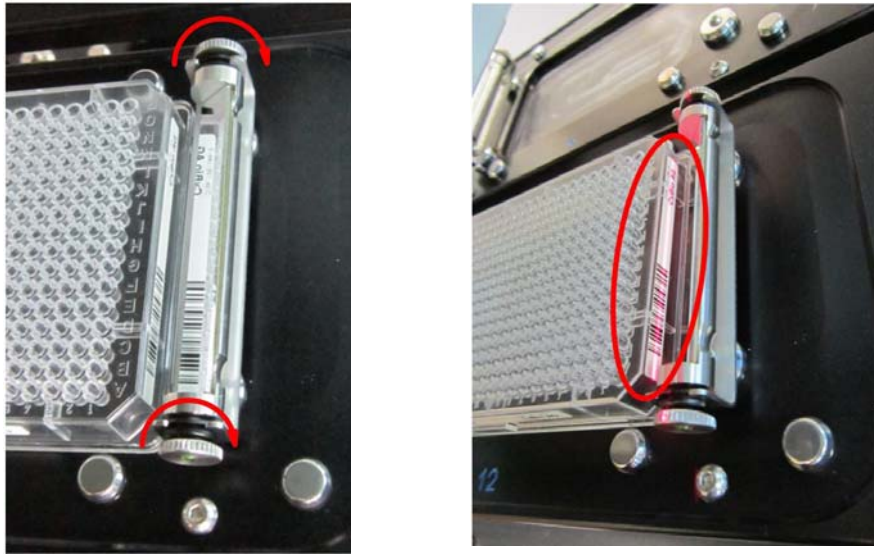**Fig. 3: Left: Adjustment of mirror
Right: Barcode reading line adjusted to middle of barcode**

Verify if red reading line crosses the microplate's barcode label horizontally. If not, loosen the two mounting screws directly at the barcode reader using Allen key 3.0 mm and tilt barcode reader horizontally. Fasten the mounting screws when completed.

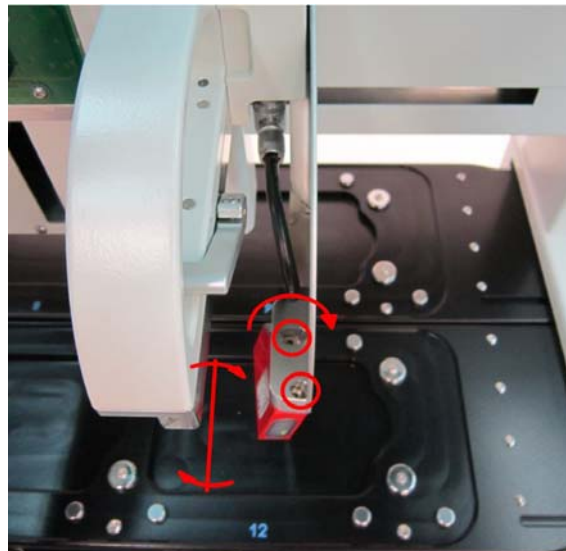
Fig. 4: Barcode reader horizontal alignment

The barcode is successfully being read if the barcode reader's decode LED (upper LED) flashes from orange to green for 200 ms, then off, then orange again repeatedly.

Close blind. Click OK on the message box. The barcode is now displayed in a new message box. Click OK to finish.

Test Barcode Reading

Put microplates with barcode labels on all deck positions to be tested. Power up instrument with pipetting head inserted but nothing else mounted to the pipetting head.

NOTE

Barcode reading requires any tips, liquid handling adapters, CyBi®-RoboTipTrays or cover magazine not being mounted to the head. The barcode reading library scripts will verify this and give an error if anything is mounted.

Start CyBio® Composer. Run empty method to calibrate drives. Add script "Test FeliX Barcodes" from the library service "CyBi®-FeliX barcode reader" to a new method. Start method. Select a data path where the script will save the barcodes to. Then select deck positions for testing barcode reading. Set the repeat counter to >1 if repeated testing is desired. The instrument will move the barcode reader to each of the selected deck position and read the barcode. A message box will show the reading result. Click OK to continue. If reading fails, you may click "Ignore" to test the next deck position. Otherwise the method will stop with an error. Check barcode quality and repeat mirror adjustment for those deck positions that failed reading the barcode.

Operation

Add script "Read barcode" from the library service "CyBi®-FeliX barcode reader" to your method at the method step where a barcode is to be obtained from the labware. Make sure any tips, liquid handling adapters, CyBi®-RoboTipTrays or cover magazine are removed before the read barcode step.

NOTE

Barcode reading requires any tips, liquid handling adapters, CyBi®-RoboTipTrays or cover magazine not being mounted to the head. The barcode reading library scripts will verify this and give an error if anything is mounted.

The library script "Read barcode" provides following parameters.

Position:

Passes deck position from main routine to Read barcode sub routine.

Barcode:

Passes read barcode from Read barcode sub routine to main routine.

Barcode time:

Passes time stamp when barcode is read from Read barcode sub routine to main routine.

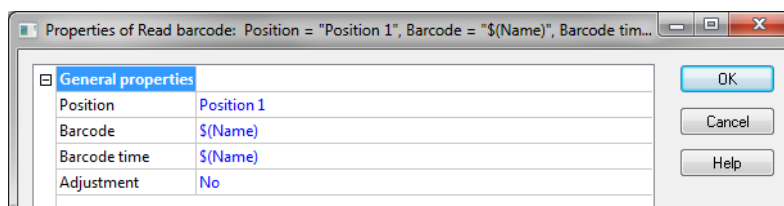
Fig. 5: Properties of Read barcode

Define appropriate variables in the main routine and assign them in the "Read barcode" sub routine parameter list.

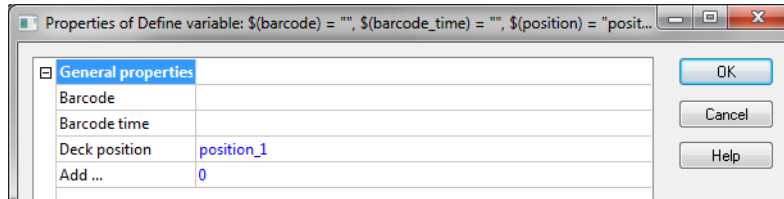
Fig. 6: Properties of Define variable

Use text as variable type for "Barcode", "Barcode time" and "Deck position"

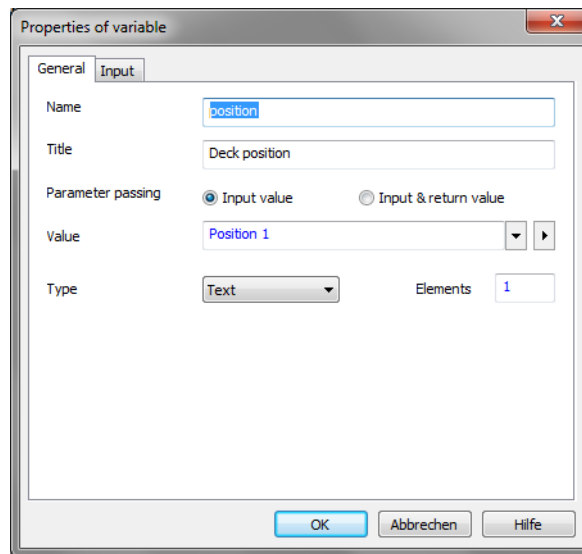
Fig. 7: Properties of variable for deck position

For "Deck position", use "Position list" as input control and select "Service": "CyBi®-Felix" and "Drive": "Virtual 3D drive".

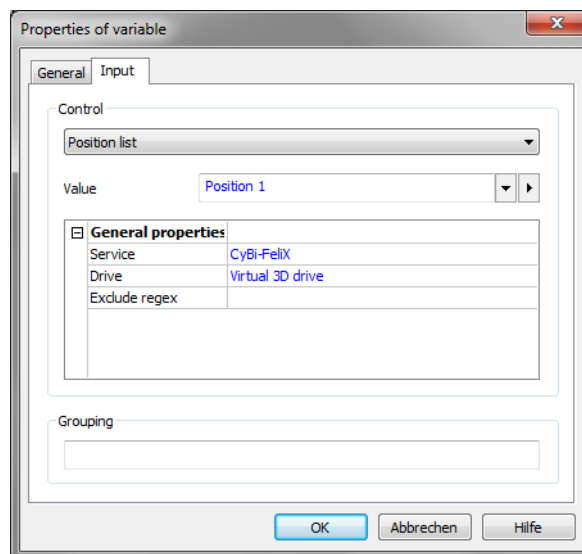
Fig. 8: Properties of variable for deck position

Assign the defined variables to the parameters of the library script "Read barcode".

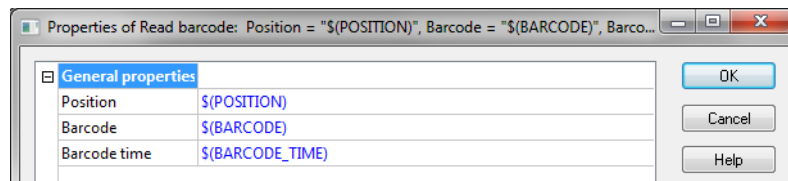
Fig. 9: Properties of Assign variable

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